


Protect every life story™ in your lab.

LaserTrack™ PH1 Cassette Printer

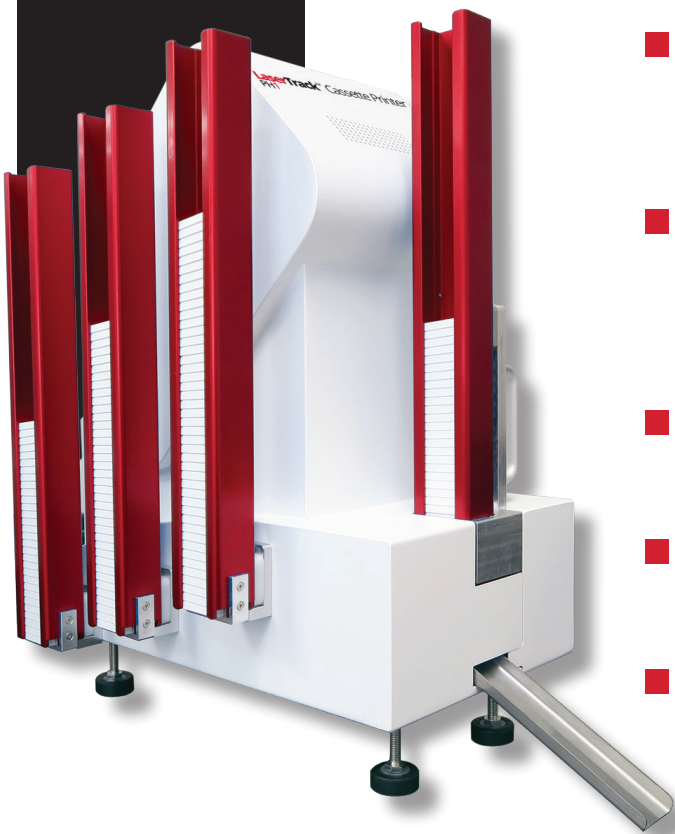

■ **Reliable, High Resolution Text & Barcodes**

Cassettes are printed using an innovative laser technology that creates the precise, clean edges and clear lines a barcode requires to produce a nearly perfect scan rate throughout the tracking process.

■ **Optimized Cassettes**

LaserTrack cassettes are specifically designed and optimized for use in LaserTrack cassette printers. This design results in effective movement within the instrument that reduces jamming and produces optimal scan rates of nearly 100%.

■ **Unsurpassed Print Speed**

The PH1 prints cassettes at a rate of one every three seconds. Four easily interchangeable magazines hold up to 80 cassettes each.

■ **Ancillary Consumables Eliminated**

The Phaser cassette is all that is required. Hot foil tape, light bulbs, ink cartridges, printheads, and filters are not necessary.

■ **Compact Footprint**

At only 15" x 8" x 18", this unit fits easily into any size work space.

Product Specifications:

Dimensions: LxDxH (15 x 8 x 18) (38.1 x 20.32 x 45.72) (in) (mm)

Weight: lbs, kg 39 lbs; 17.7 Kg

Voltage: 100 – 240V, 50/60Hz, 3A
The main power supply voltage fluctuations should not exceed ±10%.

No. of Magazines 1 (with 3 additional magazines)

Capable of interfacing with:

- Cerner
- Epic
- Sunquest
- MEDITECH
- Other LIS/LIMS systems

Characters Alphanumeric, linear barcode or data matrix ASCII codes

Network Connection USB

Electrical Conformity 2004/108/EEC 2006/95/EEC 2006/42/EEC
EM Compatibility Low Tension Devices

G1-205-CLD01 LaserTrack™ PH1 Laser Cassette Printer

Includes:

- Power Cord
- USB Cable
- User Manual
- 1 magazine with 3 spare magazines
- Operating Software

G1-2EE-9844 LaserTrack™ Collection Tray



Cassettes*

Tissue Cassettes

G1-DTP-100(color)	(WE)	(BE)	(GN)	(GY)	(PK)	(OE)	(YW)	(LC)	(BG)	(AQ)	(RE)	(OC)
-------------------	------	------	------	------	------	------	------	------	------	------	------	------

Biopsy Cassettes

G1-DBP-100(color)	(WE)	(BE)	(GN)	(GY)	(PK)	(OE)	(YW)	(LC)	(BG)	(AQ)	(RE)	(OC)
-------------------	------	------	------	------	------	------	------	------	------	------	------	------

Ordering Example: to order Green Tissue Cassettes, the part number would be G1-DTP-100GN

More Information

844.643.1129

www.general-data.com/hc

General Data Healthcare helps labs enhance productivity and improve workflow with innovative products, solutions and support that focus on delivering better patient safety & care.

©2020 General Data Healthcare, Inc. All rights reserved. Specifications are current at the time of publication; however, they are subject to change without notification. All trademarks and registered trademarks are property of General Data Healthcare, Inc., Cincinnati, Ohio, USA.



General Data Healthcare, Inc.
4043 McMann Road, Cincinnati, OH 45245
844.643.1129

Email: talktous@general-data.com
Web: www.general-data.com/healthcare

Service:
PH: 800.245.5029
Email: customersupport@general-data.com