



Protect every life story™ in your lab.

SHURTrim™ Paraffin Block Dewaxer



- ✔ **No Damage to Barcodes & Printed Text**
SHURTrim™ is the safe and practical way to remove excess paraffin wax from embedded tissue cassettes without risk of damage to barcodes and printed text.
- ✔ **Ergonomic Design Is Appealing to Both Left & Right Handed Users**
SHURTrim™ is a compact tabletop design that fits perfectly and easily in any lab. The ergonomic sloped top surface provides optimum comfort to both left and right handed users.
- ✔ **Heated Grooved Surface Heats Quickly For Fast Startup**
SHURTrim™ grooved black top surface heats 0-75°C and melts excess paraffin from embedded tissue cassettes. Large work surface allow users to de-wax multiple blocks at the same time.
- ✔ **Easier Insertion of Tissue Blocks into Microtomes**
Removal of excess paraffin wax around the cassette allows the cassette to be easily inserted into microtome quick-release clamps.
- ✔ **Drainage Cup Collects Excess Paraffin**
Melted paraffin is shaved away; and the runoff wax is collected in a disposable receptacle cup below the unit.

Product Specifications:

Dimensions: WxDxH (in)(cm)	(6 x9x7) (15 x23x18) (in)(cm)
Weight :lbs(kg)	5 (3)
Voltage:	100-120 VAC, 50/60Hz, 2A and 220-240VAC, 50/60Hz, 1A

How To Order:

ST-120	SHURTrim™ Paraffin block dewaxer with drain cup holder. Includes: 5/ea drain cups. Specifications: 100-120VAC, 50/60Hz, 2A.
ST-220	SHURTrim™ Paraffin block dewaxer with drain cup holder. Includes: 5/ea drain cups. Specifications: 220-240VAC, 50/60Hz, 1A.
ST-PCH	Paraffin Drain Cup Holder


More Information

844.643.1129

www.general-data.com/hc

General Data Healthcare helps labs enhance productivity and improve workflow with innovative products, solutions and support that focus on delivering better patient safety & care.

©2016 General Data Healthcare, Inc. All rights reserved. Specifications are current at the time of publication; however, they are subject to change without notification. All trademarks and registered trademarks are property of General Data Healthcare, Inc., Cincinnati, Ohio, USA.

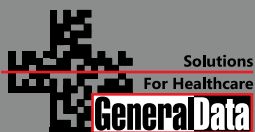

General Data Healthcare, Inc.
4043 McMann Road, Cincinnati, OH 45245
844.643.1129

Email: talktous@general-data.com
Web: www.general-data.com/healthcare

Service:
PH: 800.245.5029
Email: customersupport@general-data.com